

First edition
2017-02

AMENDMENT 1
2020-03

**Microbiology of the food chain —
Horizontal method for the detection,
enumeration and serotyping of
Salmonella —**

**Part 1:
Detection of *Salmonella* spp.**

AMENDMENT 1: Broader range of
incubation temperatures, amendment to
the status of Annex D, and correction of
the composition of MSRV and SC

*Microbiologie de la chaîne alimentaire — Méthode horizontale pour
la recherche, le dénombrement et le sérotypage des Salmonella —*

Partie 1: Recherche des Salmonella spp.

*AMENDEMENT 1: Extension de la plage de températures pour
l'incubation, amendement du statut de l'Annexe D, et correction de
la composition des milieux MSRV et SC*



Reference number
ISO 6579-1:2017/Amd.1:2020(E)

© ISO 2020



COPYRIGHT PROTECTED DOCUMENT

© ISO 2020

All rights reserved. Unless otherwise specified, or required in the context of its implementation, no part of this publication may be reproduced or utilized otherwise in any form or by any means, electronic or mechanical, including photocopying, or posting on the internet or an intranet, without prior written permission. Permission can be requested from either ISO at the address below or ISO's member body in the country of the requester.

ISO copyright office
CP 401 • Ch. de Blandonnet 8
CH-1214 Vernier, Geneva
Phone: +41 22 749 01 11
Fax: +41 22 749 09 47
Email: copyright@iso.org
Website: www.iso.org

Published in Switzerland

Foreword

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

The procedures used to develop this document and those intended for its further maintenance are described in the ISO/IEC Directives, Part 1. In particular, the different approval criteria needed for the different types of ISO documents should be noted. This document was drafted in accordance with the editorial rules of the ISO/IEC Directives, Part 2 (see www.iso.org/directives).

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights. Details of any patent rights identified during the development of the document will be in the Introduction and/or on the ISO list of patent declarations received (see www.iso.org/patents).

Any trade name used in this document is information given for the convenience of users and does not constitute an endorsement.

For an explanation of the voluntary nature of standards, the meaning of ISO specific terms and expressions related to conformity assessment, as well as information about ISO's adherence to the World Trade Organization (WTO) principles in the Technical Barriers to Trade (TBT) see www.iso.org/iso/foreword.html.

This document was prepared by Technical Committee ISO/TC 34, *Food products*, Subcommittee SC 9, *Microbiology*, in collaboration with the European Committee for Standardization (CEN) Technical Committee CEN/TC 463, *Microbiology of the food chain*, in accordance with the Agreement on technical cooperation between ISO and CEN (Vienna Agreement).

A list of all parts in the ISO 6579 series can be found on the ISO website.

Any feedback or questions on this document should be directed to the user's national standards body. A complete listing of these bodies can be found at www.iso.org/members.html.

STANDARDSISO.COM : Click to view the full PDF of ISO 6579-1:2017/Amd 1:2020

Microbiology of the food chain — Horizontal method for the detection, enumeration and serotyping of *Salmonella* —

Part 1: Detection of *Salmonella* spp.

AMENDMENT 1: Broader range of incubation temperatures, amendment to the status of Annex D, and correction of the composition of MSRV and SC

Foreword, fifth bullet

Replace the bullet with the following:

- The temperature range for incubation of non-selective and selective media has been extended from $37\text{ °C} \pm 1\text{ °C}$ to 34 °C to 38 °C without further tolerance.

4.3, second paragraph

Replace the paragraph with the following:

The RVS broth or the MSRV agar is incubated at $41,5\text{ °C}$ for 24 h and the MKTTn broth between 34 °C and 38 °C for 24 h.

4.4, last paragraph

Replace the paragraph with the following.

The XLD agar is incubated between 34 °C and 38 °C and examined after 24 h. The second selective agar is incubated according to the manufacturer's instructions.

6.3

Replace the text with the following:

6.3 Incubator, capable of operating in the range 34 °C to 38 °C .

NOTE The range 34 °C to 38 °C for incubation of media includes the use of incubators set at $35\text{ °C} \pm 1\text{ °C}$, $36\text{ °C} \pm 2\text{ °C}$ or $37\text{ °C} \pm 1\text{ °C}$.

6.6

Replace the text with the following:

6.6 Water bath, capable of operating in the range 34 °C to 38 °C.

NOTE The range 34 °C to 38 °C for incubation of media includes the use of water baths set at 35 °C ± 1 °C, 36 °C ± 2 °C or 37 °C ± 1 °C.

9.3.2, fifth paragraph

Replace the paragraph with the following:

Incubate the inoculated MKTTn broth between 34 °C and 38 °C (6.3) for 24 h ± 3 h.

9.4.2, fifth paragraph

Replace the paragraph with the following:

Incubate the XLD plates inverted between 34 °C and 38 °C (6.3) for 24 h ± 3 h.

9.4.3, second paragraph

Replace the paragraph with the following:

Incubate the XLD plates inverted between 34 °C and 38 °C (6.3) for 24 h ± 3 h.

9.5.3.2, first paragraph

Replace the paragraph with the following:

Streak the agar slant surface and stab the butt. Incubate between 34 °C and 38 °C (6.3) for 24 h ± 3 h.

9.5.3.3, first paragraph

Replace the paragraph with the following:

Streak the agar slant surface. Incubate between 34 °C and 38 °C (6.3) for up to 24 h.

9.5.3.4, first paragraph

Replace the paragraph with the following:

Inoculate just below the surface of the liquid medium. Incubate between 34 °C and 38 °C (6.3) for 24 h ± 3 h.

9.5.3.5, first line of fourth paragraph

Replace the first line of the fourth paragraph with the following:

Add one drop of toluene and shake the tube. Place the tube in a water bath and incubate between 34 °C and 38 °C (6.6) for several minutes (approximately 5 min).

9.5.3.5, fifth paragraph

Replace the paragraph with the following:

Replace the tube in the water bath and incubate between 34 °C and 38 °C (6.6) for up to 24 h.

9.5.3.6, third paragraph


Replace the paragraph with the following:

Incubate between 34 °C and 38 °C (6.3) for 24 h $\pm$ 3 h. After incubation, add 1 ml of the Kovacs reagent (B.12.2).


STANDARDSISO.COM : Click to view the full PDF of ISO 6579-1:2017/Amd.1:2020

Annex A

Replace Figure A.1 with the following figure:

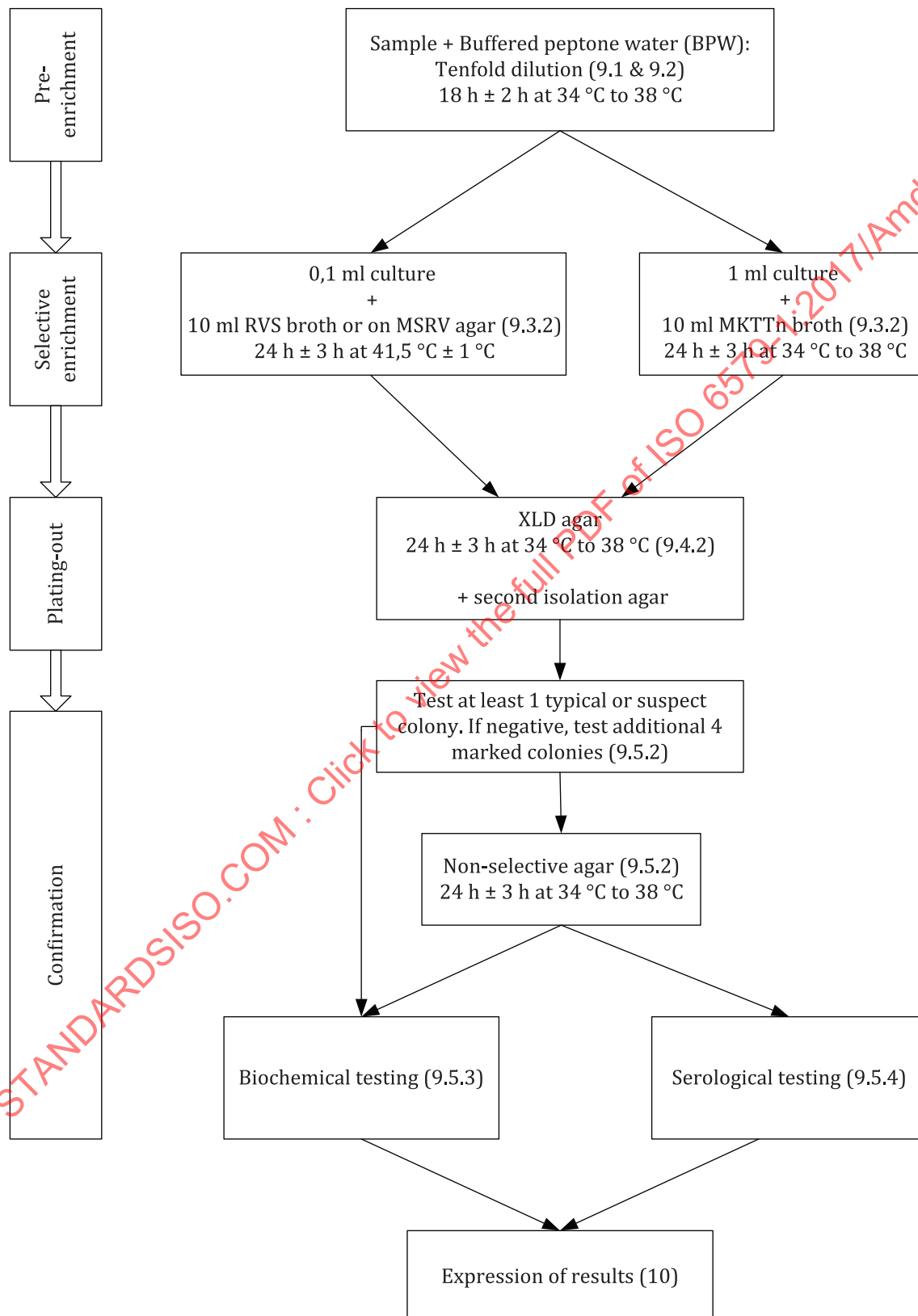

Replace Figure A.2 with the following figure:

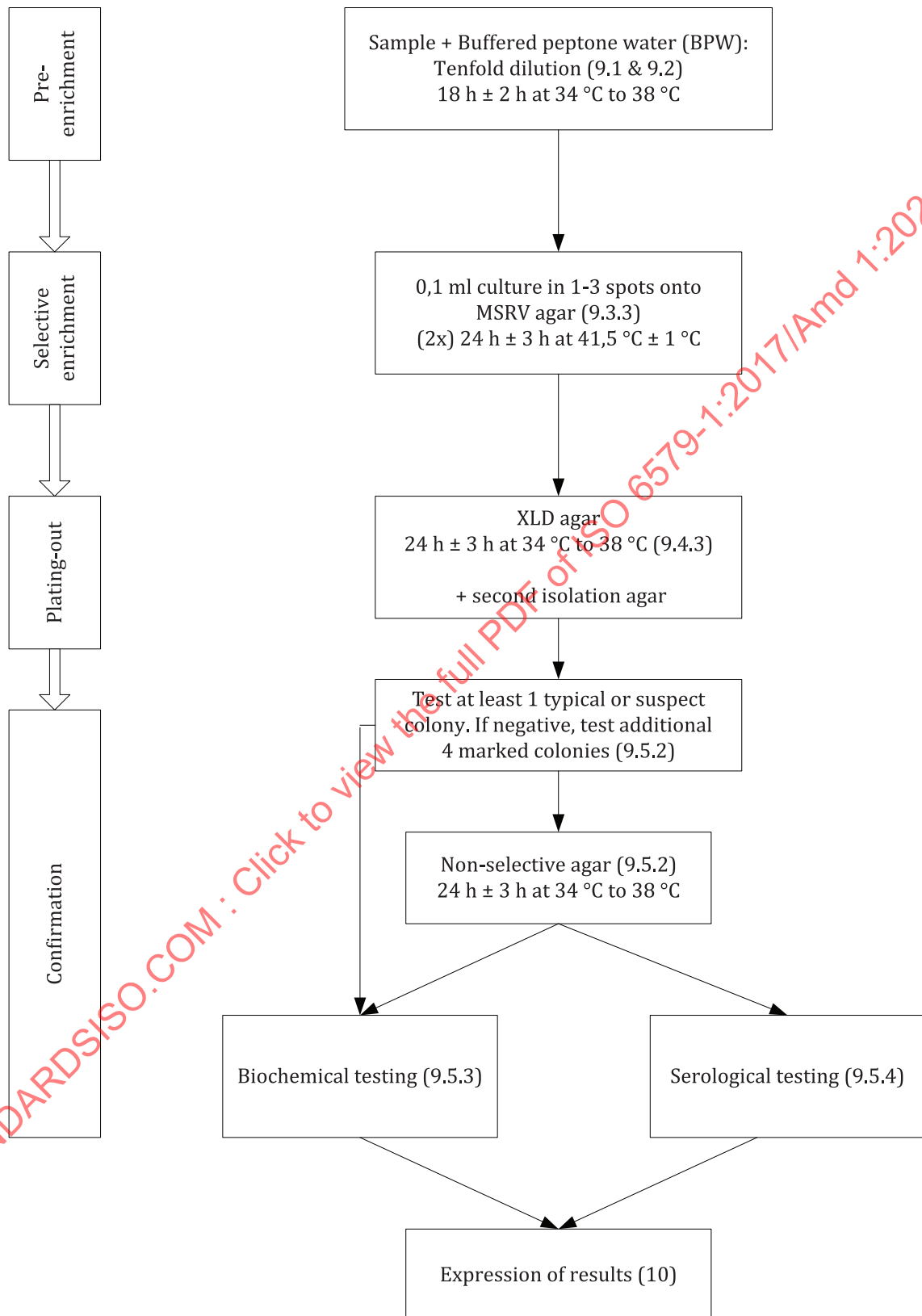

Annex B, B.4

Replace the text with the following:

B.4 Modified semi-solid Rappaport-Vassiliadis (MSRV) agar

NOTE See Reference [12].

B.4.1 Solution A

B.4.1.1 Composition

Enzymatic digest of animal and plant tissue	5,0 g
Acid hydrolysate of casein	5,0 g
Sodium chloride	8,0 g
Potassium dihydrogen phosphate (KH_2PO_4)	1,6 g
Water	1 000 ml

B.4.1.2 Preparation

Dissolve the components in the water by heating to about 70 °C, if necessary.

The solution shall be prepared on the day of preparation of the complete MSRV agar.

B.4.2 Solution B

B.4.2.1 Composition

Magnesium chloride hexahydrate ($\text{MgCl}_2 \cdot 6\text{H}_2\text{O}$)	400,0 g
Water	1 000 ml

B.4.2.2 Preparation

Dissolve the magnesium chloride in the water.

As this salt is very hygroscopic, it is advisable to dissolve the entire contents of $\text{MgCl}_2 \cdot 6\text{H}_2\text{O}$ from a newly opened container according to the formula. For instance, 250 g of $\text{MgCl}_2 \cdot 6\text{H}_2\text{O}$ is added to 625 ml of water giving a solution of total volume of 788 ml and a mass concentration of about 31,7 g per 100 ml of $\text{MgCl}_2 \cdot 6\text{H}_2\text{O}$.

The solution may be kept in a dark glass bottle with tight stopper at room temperature for at least two years.

B.4.3 Solution C

B.4.3.1 Composition

Malachite green oxalate	0,4 g
Water	100 ml

B.4.3.2 Preparation

Dissolve the malachite green oxalate in the water.

The solution may be kept in a dark glass bottle at room temperature for at least eight months.

B.4.4 Base medium**B.4.4.1 Composition**

Solution A (B.4.1)	1 000 ml
Solution B (B.4.2)	80 ml
Solution C (B.4.3)	10 ml
Agar ^a	3,0 g
^a It may be necessary to determine experimentally the concentration of agar needed for the optimal swarming of <i>Salmonella</i> (e.g. when using a batch of agar with unknown gel strength).	

B.4.4.2 Preparation

Add to 1 000 ml of solution A, 80 ml of solution B, and 10 ml of solution C and mix by agitation.

Add and suspend the agar.

Adjust the pH, if necessary, so that after sterilization, it is 5,2 (5,1 to 5,4) at 20 °C to 25 °C.

Heat to boiling with agitation. **Do not autoclave.**

Do not hold the medium at high temperatures longer than necessary.

Cool the medium to 47 °C to 50 °C (6.5).

B.4.5 Novobiocin solution**B.4.5.1 Composition**

Novobiocin sodium salt	0,05 g
Water	10 ml

B.4.5.2 Preparation

Dissolve the novobiocin sodium salt in the water.

Sterilize by filtration through a filter with a pore size of 0,2 µm.

The solution may be stored for up to four weeks at 5 °C (6.8) or in small portions (e.g. of 2 ml) at -20 °C (6.9) for up to one year.

B.4.6 Complete medium**B.4.6.1 Composition**

Base medium (B.4.4)	1 090 ml
Novobiocin solution (B.4.5)	2,2 ml

B.4.6.2 Preparation

Aseptically, add 2,2 ml of the novobiocin solution (B.4.5) to 1 090 ml of base medium (B.4.4) at 47 °C to 50 °C (6.5). Mix carefully. The final concentration of novobiocin in the complete medium is 10 mg/l.

The final pH shall be 5,2 (5,1 to 5,4) at 20 °C to 25 °C.

Pour the medium into sterile Petri dishes (6.14) up to a volume of 15 ml to 20 ml in dishes with a diameter of 90 mm.

Allow the medium to solidify before moving and handle with care.

Store the plates, **with surface upwards**, and protected from drying for up to two weeks at 5 °C (6.8) in the dark.

Do not invert the plates as the semi-solid agar is too liquid to do so.

Any plates in which the semi-solid agar has liquefied or fragmented shall not be used.

Immediately, before use and only if visible moisture is apparent, dry the surface of the agar plates carefully, for example, by placing them with the lids off and the agar surface **upwards** in a laminar air flow cabinet. Take care not to overdry the medium.

NOTE 1 The composition of MSRV agar, as described in Reference [12], contains 20 mg/l of novobiocin. However, from a scientific point of view, 10 mg/l novobiocin is preferred. Studies have shown larger migration zones on MSRV agar with a lower concentration of novobiocin^[23] and the (negative) influence of novobiocin on bacterial motility^[22].

NOTE 2 The final medium composition is enzymatic digest of animal and plant tissue 4,6 g/l, acid hydrolysate of casein 4,6 g/l, sodium chloride (NaCl) 7,3 g/l, potassium dihydrogen phosphate (KH₂PO₄) 1,5 g/l, anhydrous magnesium chloride (MgCl₂) 10,9 g/l, or magnesium chloride hexahydrate (MgCl₂·6H₂O) 23,3 g/l, malachite green oxalate 0,04 g/l, novobiocin sodium salt 0,01 g/l, and agar 2,7 g/l.

B.15, Table B.1

Replace the third row of Table B.1 with the following:

MKTn broth	Productivity	24 h ± 3 h/ 34 °C to 38 °C	<i>Salmonella</i> Typhimurium ^{c,d} <i>Salmonella</i> Enteritidis ^{c,d} + <i>Escherichia coli</i> ^d + <i>Pseudomonas aeruginosa</i>	00031 00030 00012 or 00013 00025	> 10 characteristic colonies on XLD agar or other medium of choice
	Selectivity		<i>Escherichia coli</i> ^d	00012 or 00013	Partial inhibition ≤ 100 colonies on TSA
			<i>Enterococcus faecalis</i> ^d	00009 or 00087	< 10 colonies on TSA

B.15, Table B.1

Replace the sixth row of Table B.1 with the following:

XLD agar	Productivity	24 h ± 3 h/ 34 °C to 38 °C	<i>Salmonella</i> Typhimurium ^{c,d} <i>Salmonella</i> Enteritidis ^{c,d}	00031 00030	Good growth (2) of colonies with black centre and a lightly transparent zone of reddish colour due to the colour change of the medium
	Selectivity		<i>Escherichia coli</i> ^d	00012 or 00013	Growth or partial inhibition (0–1) of yellow colonies
			<i>Enterococcus faecalis</i> ^d	00009 or 00087	Total inhibition (0)